


Rubber Fab's

BIOLOGICAL INDICATOR GASKET

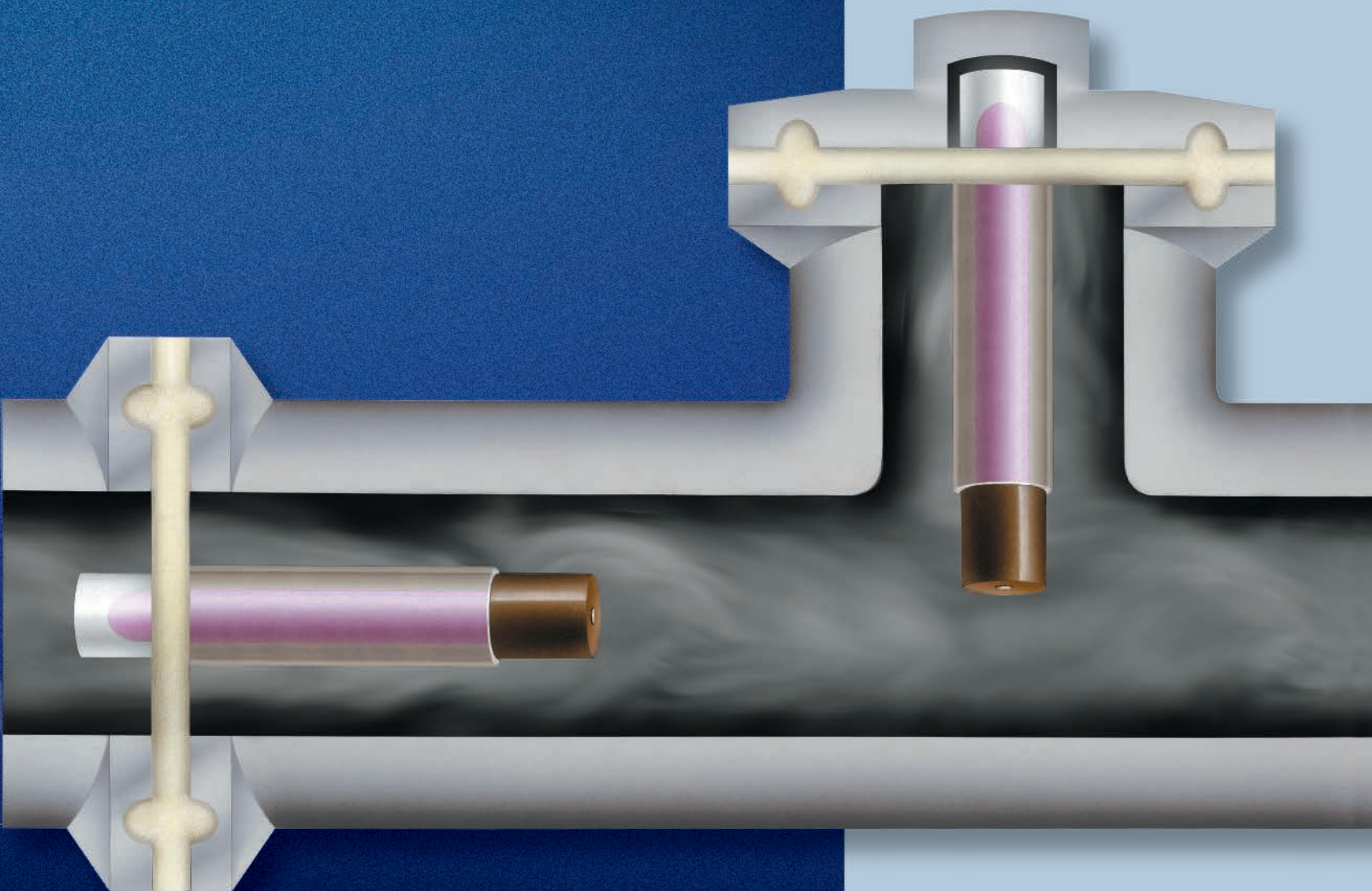

Rubber Fab

a Garlock Hygienic Technologies company

The Biological Indicator Gasket

Rubber Fab is proud to introduce a new product that enables the use of a **self contained biological indicator** in sanitary process lines. This unique Rubber Fab Smart Gasket® permits either top loading or inline positioning of an **EZTest® self contained biological indicator and/or a temperature probe**.

The great advantage of a **self contained biological indicator** over the traditional spore strip in glassine is **the time saved** when confirming spore kill. **Results are easily obtained in only 24 hours** (versus the five to seven days required with a spore strip) **and without problematic laboratory transfer**. You can feel secure knowing that the **biological indicator** will not be lost down stream in the process. The fail-safe Rubber Fab Smart Gasket holds the indicator in place for **easy retrieval**. Available for 3/4", 1.0", 1.5" and 2.0" Sanitary, Tri-Clamp process lines.

Meets High Pharmaceutical Standards Biological Indicator Gasket

- U.S. Pharmacopeia Class VI Certification
- Cytotoxicity Criteria
- CFR Title 21 Section 177.2600
- Traceability: Lot and Batch
- Certification: Lot and Batch
- ASME-BPE Standards
- USDA Standards
- 3-A Sanitary Standards
- Current Good Manufacturing Practices (CGMP)
- Metal Detectable Elastomers available
- Animal Derived Ingredient Free

EZTest®

- U.S. Pharmacopeia
- ISO 11138
- European Norm 866



Rubber Fab

a Garlock Hygienic Technologies company

26 Brookfield Drive • Sparta, NJ 07871
phone: 973-579-2959 fax: 973-579-7275

www.rubberfab.com




Additional technical and usage information about the EZTest® biological indicator and related testing products is available upon request.

Validation of Steam Systems

- Fermenters
- Pharmaceutical Tanks
- Sanitary Inline Pipe
- Pharmaceutical Hose

Smart Gasket Technology

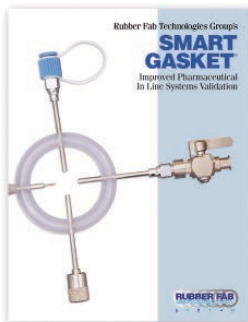

The Biological Indicator Gasket is also available with one (1) Smart Gasket® Thermocouple Port. See Rubber Fab's Smart Gasket literature on our website for more information on Rubber Fab's other inline system validation products.

© Smart Gasket and ADI Free are registered trademarks of Rubber Fab Technologies Group. Biological Indicator Gasket: United States Patent Number 6,927,058.

© EZTest is a registered trademark of SGM Biotech, Inc.

SGM BIOTECH, INC.

Rubber Fab is a member of:



Distributed by:

Rubber Fab

a Garlock Hygienic Technologies company

Biological Indicator Component

Rubber Fab is proud to introduce a technology that enables the use of a self contained biological indicator in sanitary process lines. This unique Rubber Fab Smart Gasket® permits either top loading or inline positioning of the EZTest® self contained biological indicator and/or temperature probe.



Platinum Silicone Biological Indicator Gaskets

SIZE	PART NUMBER
3/4"	42RXPX-BI-075
1"	40RXPX-BI-100
1 1/2"	40RXPX-BI-150
2"	40RXPX-BI-200

Biological Indicator Clamps (Comes with 1 Port)

SIZE	PART NUMBER
1/2" & 3/4"	CL-TH-050/075-1
1" & 1 1/2"	CL-TH-100/150-1
2"	CL-TH-200-1

Platinum Silicone Biological Indicator Gaskets w/ Smart Port

SIZE	PART NUMBER
3/4"	42RXPX-BI-075-1
1"	40RXPX-BI-100-1
1 1/2"	40RXPX-BI-150-1
2"	40RXPX-BI-200-1

Biological Indicator Short Bull Tees

SIZE	PART NUMBER
3/4"	RFB7MPS-075
1"	RFB7MPS-100
1 1/2"	RFB7MPS-150
2"	RFB7MPS-200

Biological Indicator Gasket Caps

SIZE	PART NUMBER
3/4"	BIO-CAP-075
1"	BIO-CAP-100
1 1/2"	BIO-CAP-150
2"	16AMP-200-316

